

# T-EDGE Line 10 | 11

So much More  
than just an Autoclave



- ◆ Compliance with the most strict international standards
- ◆ Clinically proven sterilization capabilities
- ◆ High- Resolution tactile touchpad
- ◆ Innovative technology
- ◆ Simple-to-use interface

# So much more than just an Autoclave

## Get your peace of mind with T-Edge

- Reliability & support - Are you looking for a safe and reliable autoclave with a long life-span? With the superior Tuttnauer quality – get your peace of mind.
- T-Edge combines durable and sturdy technology with ease-of-use.
- Our autoclaves are designed to provide high quality repeatable performance and accountability for a wide range of applications used in modern clinics of all types.
- Performance: effective cycle times that ensure the safe sterilization of all your instruments.
- Safety & traceability - Ensure highest safety with a technology that monitors all relevant process parameters. The program logs provide you certainty.
- Flexibility - With integrated water reservoir, the stand-alone concept makes the autoclave completely independent of location and easy to install.

## Optimal Chamber Capacity

Designed to use more than 75% of the chamber capacity, with 5 extra-long trays

- T-Edge 10 - 3 full (1 DIN) + 3 half (1/2 DIN) dental cassettes
- T-Edge 11 - 4 full (1 DIN) + 4 half (1/2 DIN) dental cassettes

## Compact & Elegant

- Designed for a 60 cm shelf depth
- Comfortably fits in small spaces



## Accessible & Ergonomic

- Convenient and easy loading with door opening to 110°
- Assorted water filling and drainage: top, front or automatic

T-Edge Line meets the following international standards:

### Quality Management System:

- ISO 13485
- Conforms to Medical Device Directive 93/42/EEC

### Pressure Vessel Certification route:

- ASME Code Section I & Section VIII, Div.I, PED 2014/68/EU

**Electrical Safety:** EN/IEC 61010-2-40, IEC 61010-1



## Easy Operation

- User friendly touchpad and control panel
- Automatic one-button programs, including test programs
- Control panel displays programs, cycles times, water levels, and more



## Excellent Serviceability

- Sophisticated operating system that allows fast detection of faults
- Each part can be easily removed and repaired/replaced separately

## Data & Connectivity

- Internal memory - storing up to 1,000 cycles
- Transferring Data to PC or Network via ethernet connection and USB outlet



## Ecofriendly

- Low water consumption per cycle
- Lower electricity costs with time - delay settings for machine to work off hours



**Electromagnetic compatibility:** EN 61326-1

**Verification & validation functional tests:** EN 13060, ANSI/AAMI ST55

**Software validation:** EN 62304, FDA guidance

**Usability (User experience, human**

**factor engineering) tests:** IEC 62366-1


## Documentation and Traceability



- Advanced documentation and traceability options, including barcode
- Personalized username and password



Printouts for saving hard copies (optional printer)

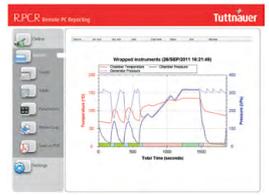

Sophisticated data analysis capabilities with R.PCR software

## Sterilisation Programs

Icon	Cycle Program	Icon	Program Name	Temp	Load Type	Type of Use
	Unwrapped 134°C		Prion 134	134°C	Solid / porous load	For storage
	Wrapped 134°C		Bowie and Dick	134°C	Chemical indicator in a product challenge device	Periodic testing as referred to in ISO 17665-1
	Unwrapped 121°C		Vacuum Test	N/A	Empty	Not applicable
	Wrapped 121°C					

## Technical Specifications

	T-EDGE10	T-EDGE11
Chamber Dim. $\phi$ x L (mm)	250 x 460	280 x 460
Chamber Volume	23 liters	27.2 liters
Dimensions W x H x D (mm)	480 x 500 x 580 mm	500 x 500 x 580 mm
Tray Dimensions (W x H x D)	421 x 189 x 20.5 mm	421 x 207 x 20.5 mm
Weight	53 kg	56 kg
Power Supply	50/60Hz, 230V-1Ph+N+GND	50/60Hz, 230V-1Ph+N+GND
Number of Trays	5	5
Connections	1 Ethernet port, 2 USB ports, WIFI, integrated automatic water filling and drain connection	
Items Supplied	Drain tube, tray holder	Drain tube, tray holder

**International Sales and Marketing**  
 Email: [info@tuttnauer.com](mailto:info@tuttnauer.com)  
[www.tuttnauer.com](http://www.tuttnauer.com)

**Tuttnauer Europe b.v.**  
 Hoeksteen 11, 4815 PR  
 PO Box 7191, 4800 GD Breda  
 The Netherlands  
 Tel: +31 765 423 510  
 Fax: +31 765 423 540  
 E-mail: [info@tuttnauer.nl](mailto:info@tuttnauer.nl)



Learn from our Experts

Join our blog [tuttnauer.com/blog](http://tuttnauer.com/blog)

